




SLM-6E Elegance DRY EYE ANALYZER

Tear film imaging unit, specially designed for dry eye examination, exquisite and precise, no loss of nasal side ring image. Large aperture ensure images high definition and color rendition. Connected to optical system with hidden wire, more controllable and stable.


Meibomian gland function evaluation

Clear observation of meibomian glands Intuitively diagnose MGD

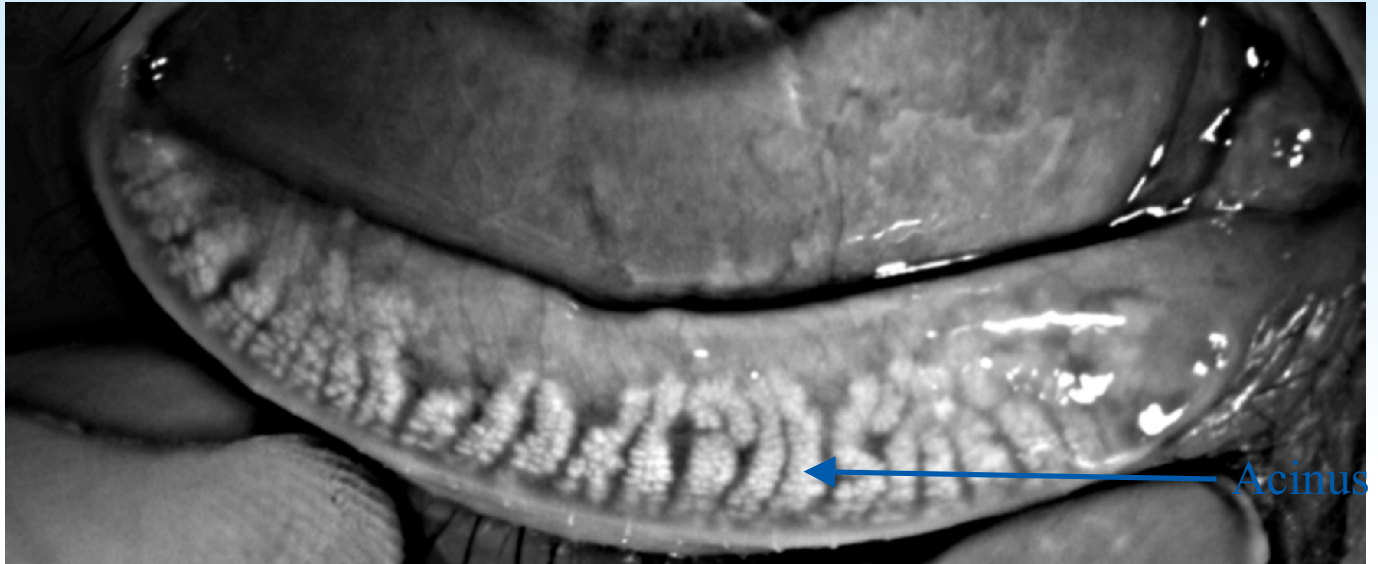

Quantitative analysis of meibomian gland deletion

With infrared imaging system and HD digital camera, the area of meibomian gland deletion can be intuitive present. Meibomian gland loss severity can be objectively evaluated by comparing with the contrast models of meibomian gland deletion level.

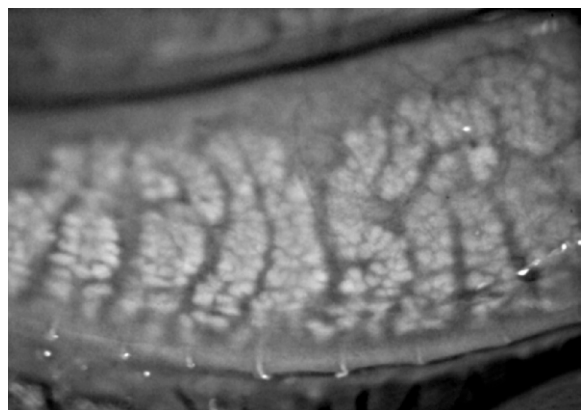

Multistep magnification observation of meibomian gland

Multistep optical magnification system support lossless image magnification. Meibomian gland distribution can be observed in the images of lower magnification, the area of missing glands can be clearly present. The detail of gland acini can be observed in the images of higher magnification, the outline is clear. It supports the clinicians to timely treat the meibomian gland lesions, diagnosis and treatment efficiency is improved.

Meibomian gland dysfunction(MGD) is the main cause of lipid layer deficiency, which leads to evaporative dry eye. It is characterized by obstruction of gland ducts, falling off the glands, or qualitative change and quantitative change of gland secretions. MGD results in tear film fast evaporation and hyperosmolarity, ocular surface inflammation, epithelium damage and discomfort.


Infrared imaging is a fast and comfortable method of meibomian gland examination, which is the only non-invasive technology of observing meibomian gland in clinic.

SLM-KD3 combines the infrared imaging technology, interference light measurement technology and multistep optical magnification technology. By taking HD images of meibomian gland and eyelid margin and dynamic video of lipid layer, the loss of meibomian gland is quantitatively classified, the eyelid margin and gland opening is intuitively observed, and the lipid layer is intuitively observed and graded. The meibomian gland function can be comprehensively analyzed accordingly, which can be the objective basis of dry eye main cause MGD diagnosis.

Eyelid margin analysis

Eyelid margin change is one of the typical characters of MGD, including eyelid margin morphologic change and meibomian gland opening change.


SLM-KD3 has multistep optical magnification system, meeting the observation requirement of eyelid margin morphologic change and meibomian gland opening change. Meibomian gland and gland opening distribution can be observed in the images of lower magnification, the details of eyelid margin and gland opening change can be observed in the images of higher magnification. The images provide the intuitive diagnosis basis for eyelid margin abnormality or MGD.

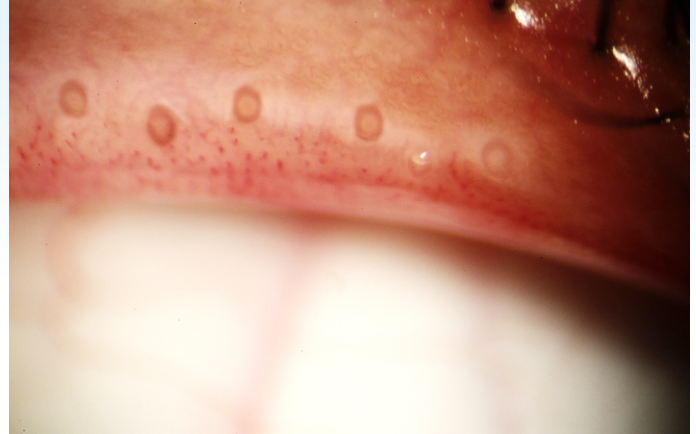

Lipid layer analysis


SLM-KD3 analyzes the thickness, quality and quantity of lipid layer by comparing with graded contrast models, diagnose MGD combining with meibomian gland and eyelid margin analysis.

Specially designed tear film imaging unit enlarges the aperture and reduces the light loss, increase the image resolution and color rendition. The lipid layer morphology and coating process can be dynamically observed by HD videos, which have exquisite frame and strong hierarchical effect.

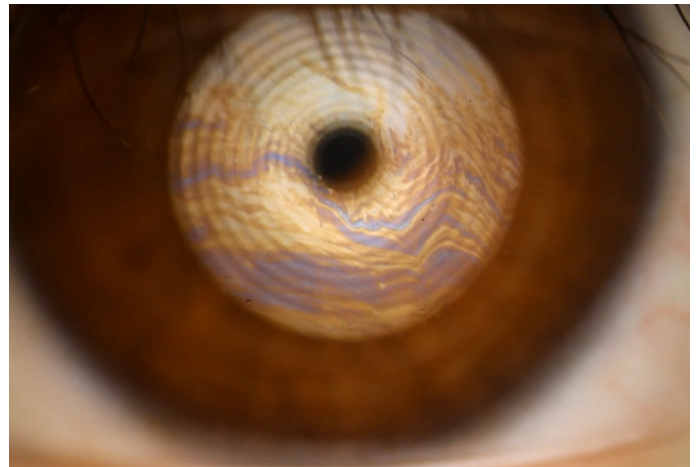

Dry eye disease comprehensive examination

All-round examination items Comprehensive dry eye evaluation

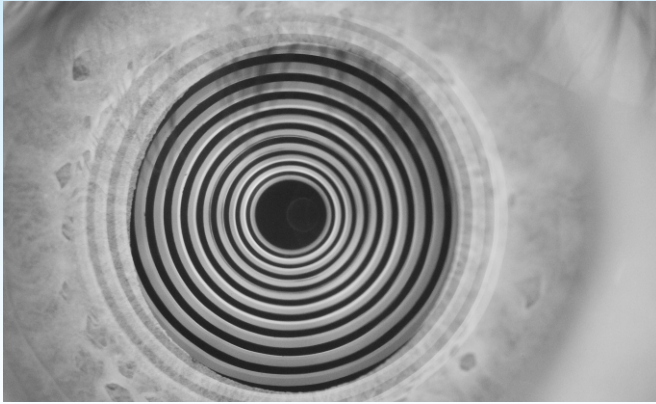


Non-invasive break up time (NIBUT)

Dry eye disease (DED) is characterized by tear film instability. SLM-KD3 quantitatively evaluates the state of tear film stability by identifying the ring image changes on the tear film.

Tear film imaging unit utilize the infrared light, which is non-invasive and non-stimulative. It directly face to the eye being examined, the ring image is complete and the nasal side ring image have no defect, which ensure the complete video record.

Compare with fluorescein staining, non-invasive examination provide more comfortable experience for the patients, automatic timing system is more accurate and efficient, and operation process is simplified.

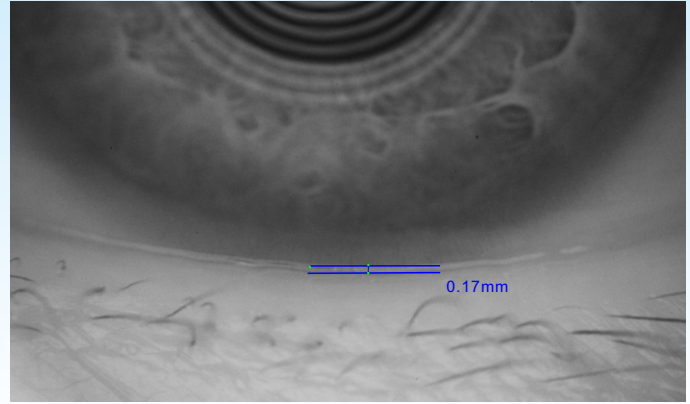


Tear meniscus height (TMH)

Compare with Schirmer test, SLM-KD3 shorten the examination time by more than 10 times, and the result is more accurate.

Tear film imaging unit utilize the infrared light, which is non-invasive and non-stimulative. Specially designed tear film imaging unit enlarges the aperture and reduces the light loss, increase the image resolution.


Multistep optical magnification system support image optical magnification for more detailed observation of tear meniscus morphology. If the tear meniscus is discontinuous or uneven, conjunctivochalasis (dry eye complication) or palpebral margin abnormal (MGD) can be simultaneously examined.

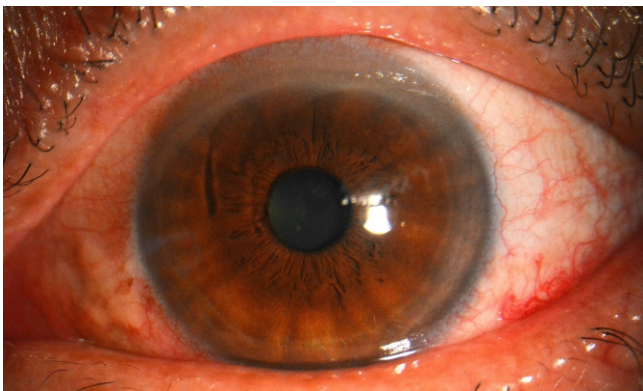

Bulbar redness scan (R-scan)

Common bulbar redness include conjunctival congestion , ciliary congestion and mixed congestion, which probably caused by dry eye disease or conjunctive. In dry eye examination, bulbar redness is preferentially considered to be aqueous deficient dry eye (ADDE).

SLM-KD3 capture images in high definition, automatically grade the conjunctival congestion and ciliary congestion and evaluate the severity of bulbar redness.

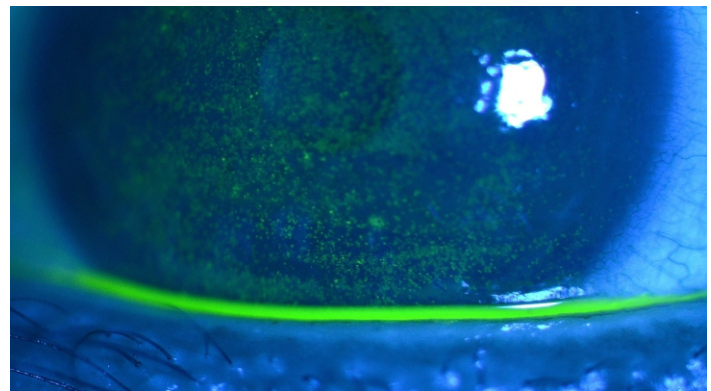



Corneal staining

Assess the corneal epithelial completeness, check mucous deficient dry eye or corneal inflammation.

Unique optical magnification system support lossless image magnification. Five blue filter combinations enhance the effect of fluorescein staining, cornea details can be observed clearly, which spares no details of small corneal lesion.

Anterior segment examination on same device


More practical function More efficient examination

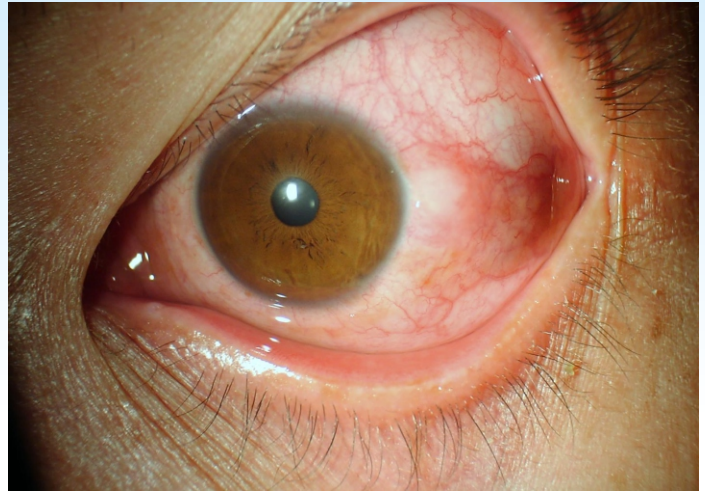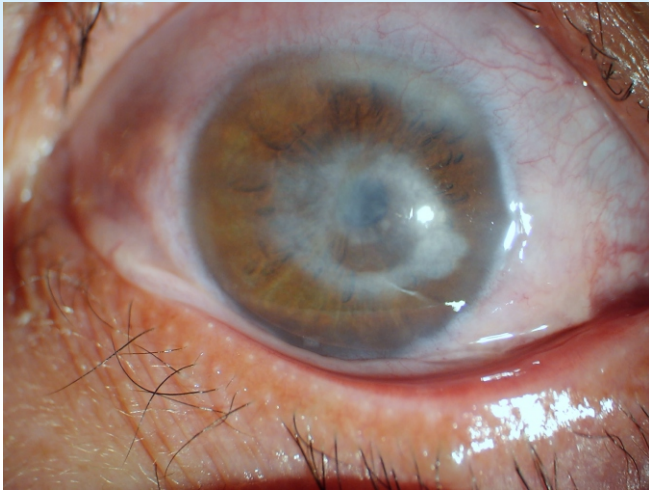

Dry eye complications examination

Without changing the device, dry eye complications can be rapidly checked, including cornea damage (cornea scarification, cornea inflammation, cornea ulcer, cornea perforation), corneal nebula, conjunctivochalasis, which improving the checking efficiency and patient's satisfaction.

Binocular eyepiece for better observation with stereoscopic view. Multistep optical magnification system support lossless image magnification. Built-in contrast-enhanced filter for more clear observation of lesion detail.

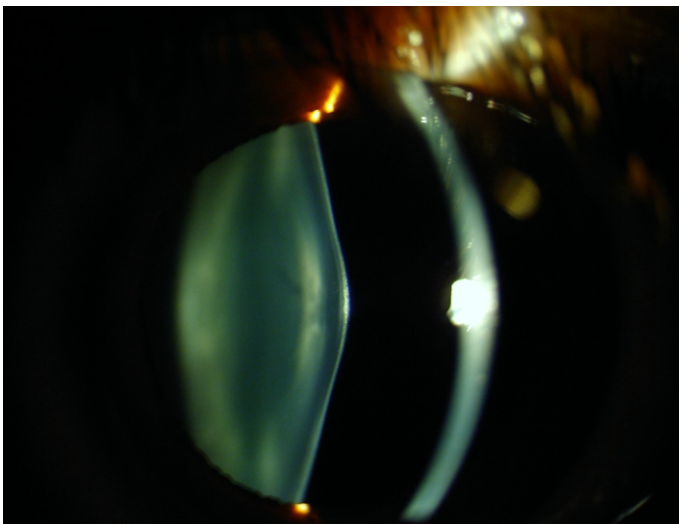

Anterior segment examination


All of anterior segment examinations can be done on same device, more efficient. Anterior segment examinations include but not limited to:

Ocular surface disease, measurement of cornea diameter and pupil diameter and lesion area, corneal lesions, lens lesions, 1/3 forepart of vitreous lesions.

Examinations of the fundus and posterior segment of vitreous can be done with other additional accessories, corneal contact lens matching examination for identifying the contact lens contraindications.


SLM-6E Elegance

Objective analysis of dry eye disease
Integrated comprehensive platform

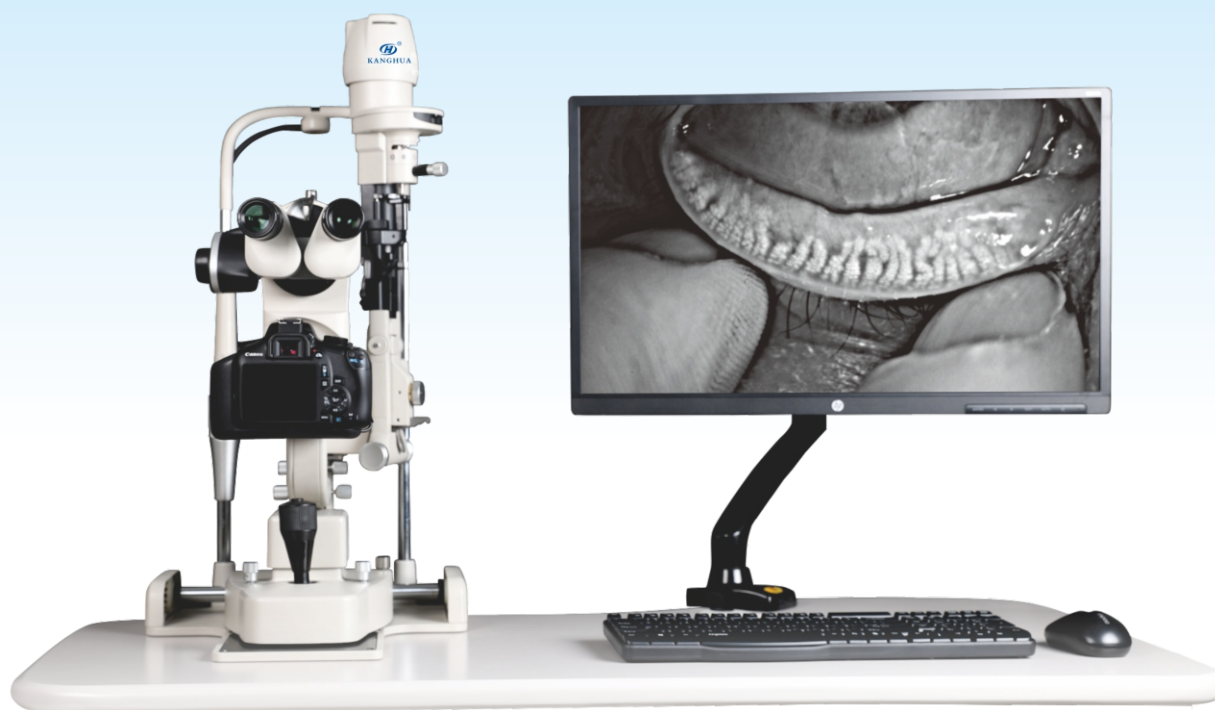

Tear film imaging unit, optical system, digital imaging system and precision mechanical platform be connected by hidden wires with special electronic technology, the device without external wire is more clean and tidy, avoid the malfunctions caused by twisted or damaged wires, improve the operation continuity and precision.

**High definition image,
Is the basis of all the digitization.**










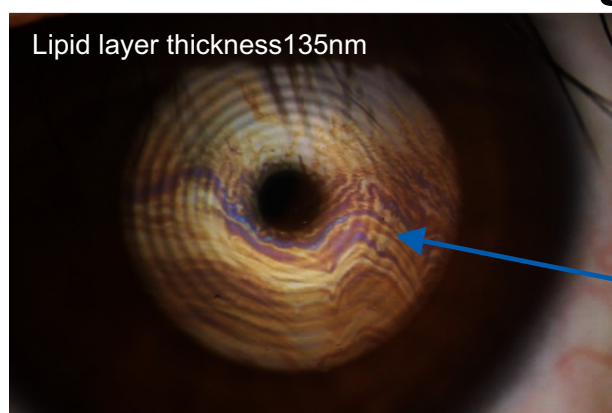


LLT(nm)	Monsell color system	color
30		White
45		Grey (White)
60		Grey
75		Grey (yellow)
90		Yellow
105		Yellow (brown)
120		Brown(yellow)
135		Brown
150		Brown(blue)
165		Blue(brown)
180		Blue

Comprehensive analysis report

Automatically generate reports Full presentation of results

Automatically generate reports* according to the examination results without manually filling in. illustrated report is easy to understand, operation process is simpler. Visible and intuitive report facilitate communication with patients.

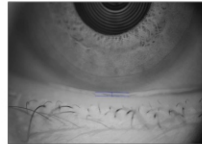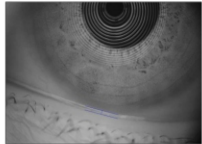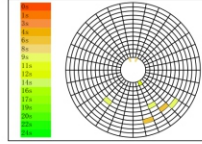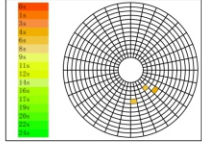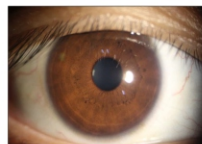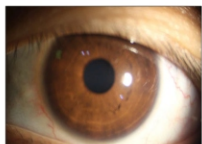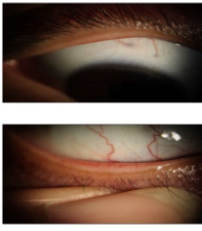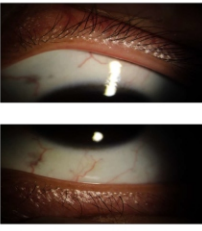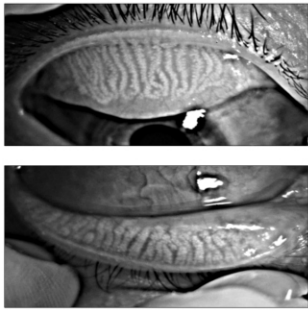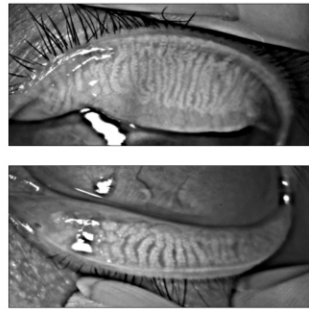Comprehensive Analysis Report			
Name: KH Sex: Male Age: 20 ExamID: 20190729004S Date: 2019/7/29 12:28:08			
	R	L	Reference value
TMH			R: 0.24mm L: 0.25mm Reference: Infrared >= 0.2mm Visible >= 0.3mm
NIBUT			R:Breakup(1st): 7.60s Breakup(avg.): 9.40s Observation time: 13.20s L:Breakup(1st): 5.70s Breakup(avg.): 6.00s Observation time: 8.90s Reference Level 0 Normal, Breakup(1st): 10 s, Breakup(avg.): 14 s. Level 1 Critical, Breakup(1st): 6-9 s, Breakup(avg.): 7-13 s. Level 2 Dry eye, Breakup(1st): 5 s, Breakup(avg.): 7 s.
Lipid layer	Level: 4	Level: 4	Reference: Level 1:<15nm Level 2:=15nm Level 3:=30nm Level 4:[30-80nm] Level 5:=80nm Level 6:[80-120nm] Level 7:[120-160nm]
R-Scan			R: Conjunctival level: 1.24 Ciliary level: 1.19 L: Conjunctival level: 1.52 Ciliary level: 1.70 Reference: <=2 Normal >2 Abnormal
Gland opening			R: Upper: Normal Nether: Normal L: Upper: Normal Nether: Normal Reference: gland openings 4 levels 1.Normal: blepharon is clear and transparent 2.Mild: gland openings bulge 3.Moderate: Eyelid margin of mucous membrane disappear and cornification 4.Severe: Eyelid irregularity,gland opening disappear,new vessels
MGD			R: Upper: 0.40 L: Upper: 0.70 Nether: 0.50 Nether: 0.40 Reference: score 0: Normal score 1: Loss <= 1/3 score 2: Loss 1/3-2/3 score 3: Loss > 2/3
Doctor: 2019-07-29 12:35			

* Comprehensive report provides all-sided results, and itemized reports present detailed images.

Dry eye	Placido ring (23 rings)	☆
	Tear meniscus height measurement	Visible light or infrared light selectable, optical magnification
	NIBUT	Visible light or infrared light selectable, duration ≤25 sec, auto calculate the first break up time and average break up time
	Lipid layer grade evaluation	Grade with contrast models, measure lipid layer thickness
	Quantitative analysis of meibomian gland deletion	Imaging with infrared, quantitative analysis of deletion, optical magnification
	Eyelid margin analysis	Gland opening grading, optical magnification
	Bulbar redness analysis	Auto grading conjunctival congestion and siliary congestion severity
Optics	Resolution ≥ 1800•N lp/mm	☆
	High sensitivity	☆
	High-point eyepiece	☆
	Magnification 6.3× 10× 16× 25× 40×	☆
Illumination	LED (white/infrared)	☆
	Slit width 0~14mm	☆
	Slit length 0.2~14mm	☆
	Slit angle 0°~180°	☆
	Inclination 5° 10° 15° 20°	☆
	Aperture diameter Φ0.2~14mm	☆
	Colorless/ Heat absorption/ Grey/ Redfree/ Cobalt blue	☆
	Five blue filter combinations	☆
	Contrast-enhanced filter	☆
Digital	1800 Megapixel DSLR camera	☆
	Auto exposure	☆
	Auto optimization	☆
	Auto eye position identification	☆
	Coaxial background light	☆
	1080P Video	☆
Others	Dry eye complications examination	☆
	Anterior segment examination	☆

Chongqing Kanghua Ruiming S&T Co.,Ltd.

Add: No.5, Road 1, Tongjiaxi Industrial Park, Beibei, Chongqing, China

Service hotline: +86-23-63227333

Zip code: 400700

E-mail: kanghuaa@kanghuaa.com

Web: www.kanghuaa.com



Company website

WE MAKE VISION MORE WONDERFUL